

ULTRA BODY

TOTAL BODY
VOLUMIZING SOLUTION

Your new experience, your last choice


HA BODY FILLER 60cc / 1vial

Natural Lines
& Volume

Voluminous
Body Lines

Cross-linked BODY Filler

ULTRA BODY is a cross-linked HA dermal filler for your body and it can restore, create volume and also fill your deepest wrinkle also.

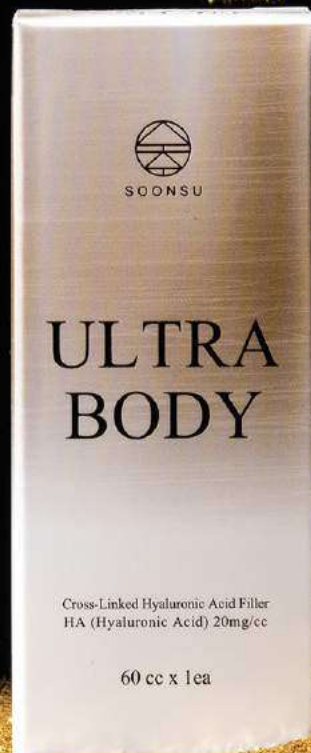

60 cc x 1ea

Cross-Linked Hyaluronic Acid Filler
HA (Hyaluronic Acid) 20mg/cc

BODY


Cross-linked BODY Filler

ULTRA BODY is a cross-linked HA dermal filler for your body and it can restore, create volume and also fill your deepest wrinkle also.



THREE-dimensional Structures

ULTRA BODY is a highly viscous filler which has three-dimensional structures and it is the result of combining HA with BDDE.

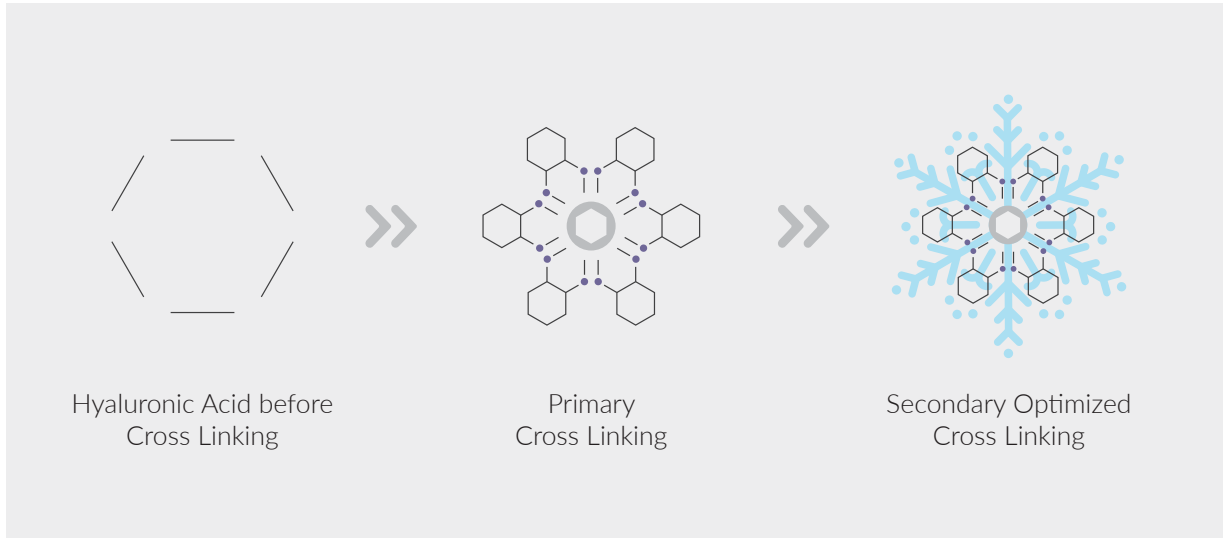



Differentiated Viscosity

ULTRA BODY keeps your shape intact.




Purity makes Safety

ULTRA BODY only uses high purity(100% Cross-linked HA) and Non-animal polymeric hyaluronic acid(EP grade). ULTRA BODY added lidocaine listed in the European EDQM and FDA. Due to 3 Times refining processing it does not leave BDDE residues . So, those main quality management processes have minimized the risk of side-effects.

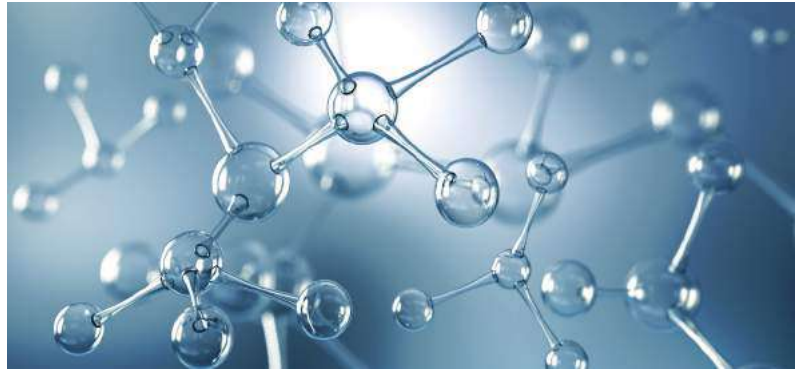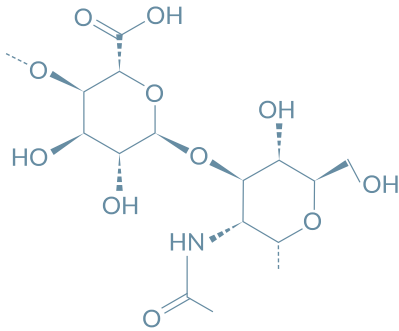


There is no regret

ULTRA BODY is effectively absorbed by injecting Hyaluronidase into the area where the patients are not satisfied.

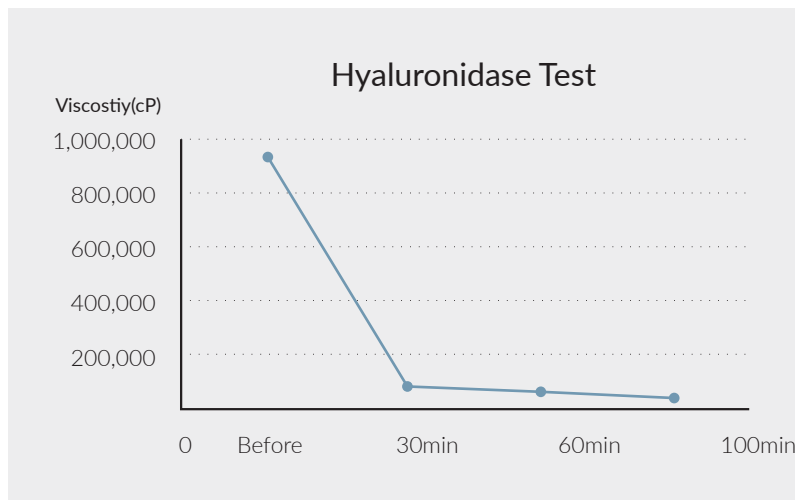

No pain ALL gain

The visual analogue score(VAS) as below.

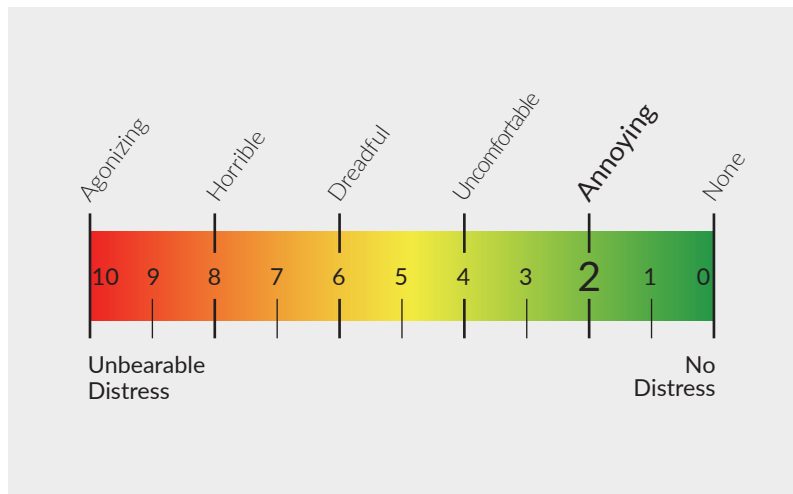

Test Result

Test Item	Standard	Result
Property and State	Colorless and clear liquid with viscosity, There must not be any foreign substances.	Conformity
Remaining BDDE	Less than or equal to 2.0 ppm	Lower than detection limits
pH	6.8 ~ 7.5	7.0
Injected Power	30N below	14.9
Capacity	More than or equal to 1.0mL	1.07 ~ 1.09 mL
Heavy Metal	Less than or equal to 10.0 ppm	Less than equal to 10 ppm
Osmotic Pressure	200 ~ 400 mmol / kg	300 mmol / kg
Asepsis / Toxicity	Less than or equal to 13.5 EU / mL	Less than or equal to 0.01 EU / mL
Aseptic / Nontoxic	Aseptic / Nontoxic	Conformity

Specification

HA Composition	20mg/cc
Volume of Vial	60cc
Inject Area	Breast, Buttocks, Local depression Area to create body volume
Lidocaine	Included
Duration	13month
Storage	2°C~30°C

How to use

USAGE 1

1. Prepare 18G, 6~10ml syringe to enter the bottle.
2. Extract the liquid and use it right away.

USAGE 2

1. Prepare 18G, 10ml x 7syringe with physiological saline.
2. Take 5ml~12ml of physiological saline and inject to the bottle.
3. After mixing, extract the liquid and use it right away.



Description

ULTRA BODY is a sterile cross-linked hyaluronic acid gel,colorless,trans parent and odorless, processed from non-cross-linked hyaluronic acid supplied by a non-animal origin.

ULTRA BODY is presented in a single-use syringe, sterilized and ready for injection.




Specification

HA Composition	20mg/cc
Volume of Vial	10ml
Inject Area	Breast, Buttocks, Local depression Area to create body volume
Lidocaine	Included
Duration	13month
Storage	2°C~30°C




Breast Before & after

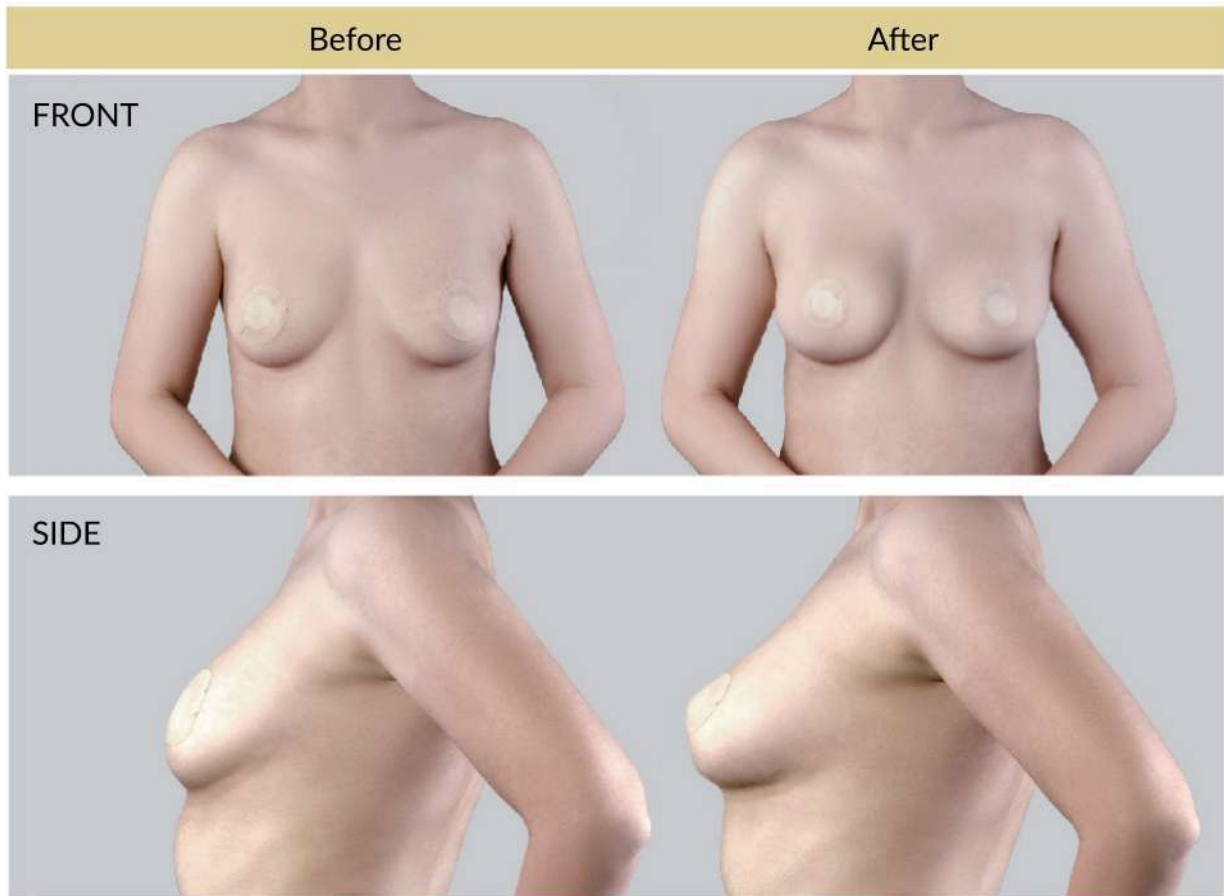


WEAR UNDERWEAR



Butt Before & after

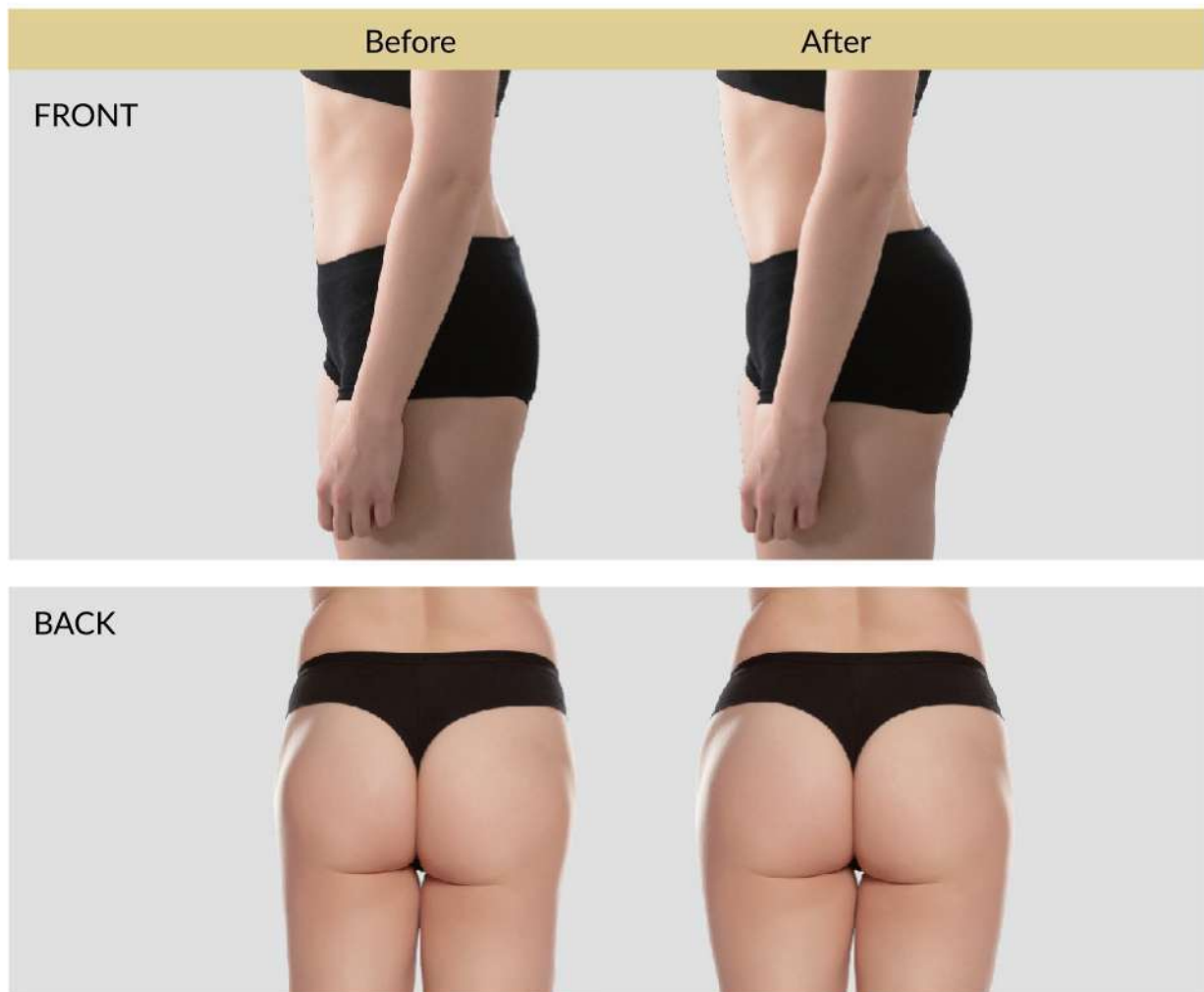

SHAPE OF BUTT

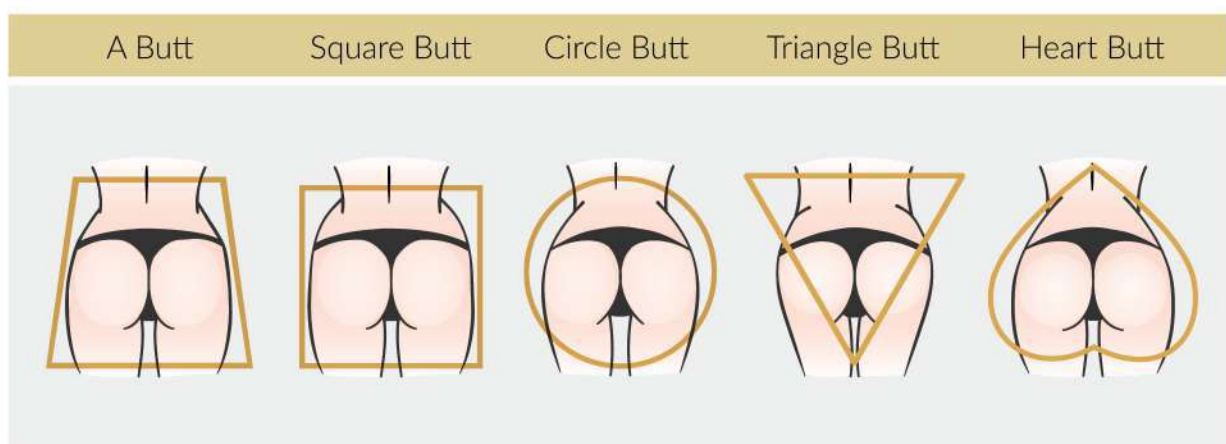



ULTRA BODY